



Non WHO Reference Material
TRANSFORMING GROWTH FACTOR BETA-1 (Bovine)
NIBSC code: 89/516
Instructions for use
(Version 6.0, Dated 23/04/2013)

This material is not for in vitro diagnostic use.

1. INTENDED USE

As a biological reference standard for transforming growth factor-beta1 (TGF-β1)

2. CAUTION

This preparation is not for administration to humans or animals in the human food chain.

The preparation contains material of human origin, and either the final product or the source materials, from which it is derived, have been tested and found negative for HBsAg, anti-HIV and HCV RNA. The bovine TGF-β1 was sourced from a country where bovine spongiform encephalopathy (BSE) has not been found. As with all materials of biological origin, this preparation should be regarded as potentially hazardous to health. It should be used and discarded according to your own laboratory's safety procedures. Such safety procedures should include the wearing of protective gloves and avoiding the generation of aerosols. Care should be exercised in opening ampoules or vials, to avoid cuts.

3. UNITAGE

The assigned potency is 1,500 arbitrary units per ampoule.

4. CONTENTS

Country of origin of biological material: United Kingdom.
Each ampoule contains the residue after freeze drying of 1ml of phosphate buffered saline that contained:

Transforming growth factor-beta1 (TGF-β1), bovine, 80ng
6mg human serum albumin


5. STORAGE

For economy of use, it is recommended that the solution be sub divided into several small aliquots and stored at -40°C or below. Unopened ampoules should be stored at -20°C.

Please note: because of the inherent stability of lyophilized material, NIBSC may ship these materials at ambient temperature.

6. DIRECTIONS FOR OPENING

Tap the ampoule gently to collect the material at the bottom (labelled) end. Ensure ampoule is scored all round at the narrow part of the neck, with a diamond or tungsten carbide tipped glass knife file or other suitable implement before attempting to open. Place the ampoule in the ampoule opener, positioning the score at position 'A'; shown in the diagram below. Surround the ampoule with cloth or layers of tissue paper. Grip the ampoule and holder in the hand and squeeze at point 'B'. The ampoule will snap open. Take care to avoid cuts and projectile glass fragments that enter eyes. Take care that no material is lost from the ampoule and that no glass falls into the ampoule.

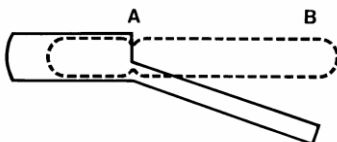

Side view of ampoule opening device containing an ampoule positioned ready to open. 'A' is the score mark and 'B' the point of applied pressure.

7. USE OF MATERIAL

No attempt should be made to weigh out any portion of the freeze-dried material prior to reconstitution

Dissolve the total contents of the ampoule with 0.5ml of sterile distilled water. Rinse the ampoule with about 0.4ml of sterile distilled water and make up the total volume to 1.0ml with sterile distilled water. This solution should now contain TGF beta 1 at a concentration of 1,500 Units/ml. Use carrier protein where extensive dilution is required.

8. STABILITY

Reference materials are held at NIBSC within assured, temperature-controlled storage facilities and they should be stored on receipt as indicated on the label. Accelerated degradation studies have indicated that this material is suitably stable, when stored at -20°C or below, for the assigned values to remain valid until the material is withdrawn or replaced. These studies have also shown that the material is suitably stable for shipment at ambient temperature without any effect on the assigned values. Once reconstituted, diluted or aliquoted, users should determine the stability of the material according to their own method of preparation, storage and use. Users who have data supporting any deterioration in the characteristics of any reference preparation are encouraged to contact NIBSC.

9. REFERENCES

This standard was produced under WHO guidelines as cited in the WHO Technical Reports Series 800, 1990, Annex 4.

10. ACKNOWLEDGEMENTS

N/A

11. FURTHER INFORMATION

Further information can be obtained as follows;
This material: enquiries@nibsc.org
WHO Biological Standards:
<http://www.who.int/biologicals/en/>
JCTLM Higher order reference materials:
<http://www.bipm.org/en/committees/jc/jctlm/>
Derivation of International Units:
http://www.nibsc.org/standardisation/international_standards.aspx
Ordering standards from NIBSC:
<http://www.nibsc.org/products/ordering.aspx>
NIBSC Terms & Conditions:
http://www.nibsc.org/terms_and_conditions.aspx

12. CUSTOMER FEEDBACK

Customers are encouraged to provide feedback on the suitability or use of the material provided or other aspects of our service. Please send any comments to enquiries@nibsc.org

13. CITATION

In all publications, including data sheets, in which this material is referenced, it is important that the preparation's title, its status, the NIBSC code number, and the name and address of NIBSC are cited and cited correctly.

14. MATERIAL SAFETY SHEET



| Physical and Chemical properties | |
|--|---|
| Classification in accordance with Directive 2000/54/EC, Regulation (EC) No 1272/2008: Not applicable or not Physical appearance: White powder | Corrosive: No |
| Stable: Yes | Oxidising: No |
| Hygroscopic: Yes | Irritant: No |
| Flammable: No | Handling: See caution, Section 2 |
| Other (specify): | Contains human and bovine source material |
| Toxicological properties | |
| Effects of inhalation: | Not established, avoid inhalation |
| Effects of ingestion: | Not established, avoid ingestion |
| Effects of skin absorption: | Not established, avoid contact with skin |
| Suggested First Aid | |
| Inhalation: | Seek medical advice |
| Ingestion: | Seek medical advice |
| Contact with eyes: | Wash with copious amounts of water. Seek medical advice |
| Contact with skin: | Wash thoroughly with water. |
| Action on Spillage and Method of Disposal | |
| Spillage of ampoule contents should be taken up with absorbent material wetted with an appropriate disinfectant. Rinse area with an appropriate disinfectant followed by water. Absorbent materials used to treat spillage should be treated as biological waste. | |

15. LIABILITY AND LOSS

In the event that this document is translated into another language, the English language version shall prevail in the event of any inconsistencies between the documents.

Unless expressly stated otherwise by NIBSC, NIBSC's Standard Terms and Conditions for the Supply of Materials (available at http://www.nibsc.org/About_Us/Terms_and_Conditions.aspx or upon request by the Recipient) ("Conditions") apply to the exclusion of all other terms and are hereby incorporated into this document by reference. The Recipient's attention is drawn in particular to the provisions of clause 11 of the Conditions.

16. INFORMATION FOR CUSTOMS USE ONLY

| |
|---|
| Country of origin for customs purposes*: United Kingdom * Defined as the country where the goods have been produced and/or sufficiently processed to be classed as originating from the country of supply, for example a change of state such as freeze-drying. |
| Net weight: 1g |
| Toxicity Statement: Toxicity not assessed |
| Veterinary certificate or other statement if applicable. |
| Attached: No |